









| | | | | |
|--|--|---|---|--|
| DLS 8 Samples 02 TR Report | | INTEGRAL VIALAB | | |
| Program Name | | DLS 8 Samples 02 TR.iaa | | |
| Program Name (on pipette) | | DLS-8_25JAN_07 | | |
| Last Saved Date: | | 25. Jan 2024 | | |
| Last Save Operator: | | NDziuba | | |
| Instrument - Serial Number | | 0020050843 | | |
| Pipette - Serial Number | | 0007018881 | | |
| Tip Type (PN 6565) Lot Nr.: | | | | |
| Run Operator: | | NDziuba | | |
| Run Date: | | 25.Jan.2024 | | |
| Run Start Time: | | 16:54 | | |
| Run End Time: | | 17:08 | | |
| Notes: | | 16:54:20 : Run started 16:54:23 : Repeat Dispense (Step 02) 16:55:58 : Repeat Dispense (Step 03) 16:57:15 : Wait for user input. Ready to continue protocol 16:57:16 : Message (Step 04) 16:57:44 : Run continued 16:57:45 : Serial Dilution (Step 05) 17:08:01 : Run finished | | |
| Signature: | | | | |
| Overview Method | | | | |
|  VOYAGER - 50µl - 8CH | <div>1 Initial Volumes</div>  | <div>2 Repeat Dispense</div>  | <div>3 Repeat Dispense</div>  | <div>4 Message</div>  |
| | <div>5 Serial Dilution</div>  | | | |
| Total Time: | | 11 min 12 sec | | |
| Total Tip Consumption: | | 24 | | |
| Deck Layout | | | | |



Pipette & Deck


| Labware | Name | Manufacturer | Part Number |
|-------------|---|--------------|-------------|
| Pipette | VOYAGER 50 µl 8 channels | INTEGRA | 4726 |
| Pipette Tip | 50/125 µl GripTip, Sterile, Filter, Low retention | INTEGRA | 6565 |
| Deck | 3 Position Universal Deck | INTEGRA | 4520 |

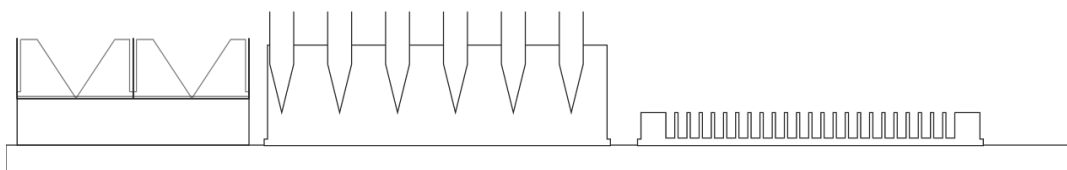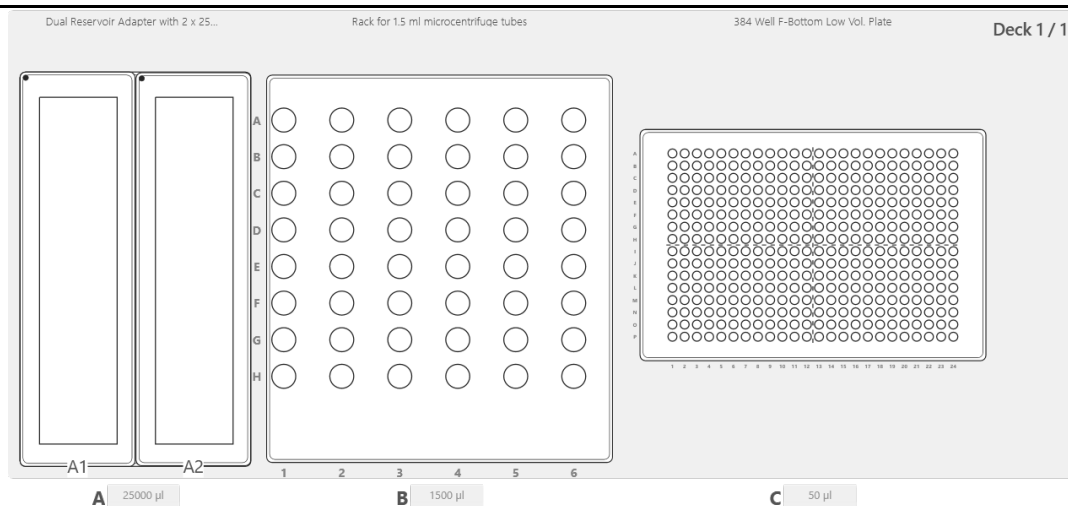Deck Labware

| Deck Position | Labware | Name | Manufacturer | Part Number | Description |
|---------------|--------------|--|--------------|--|---|
| A | COMBI System | Dual Reservoir Adapter with 2 x 25 ml Reservoirs | INTEGRA | 4547 | Dual Reservoir Adapter (PN 4547) with 2 x 25 ml Multichannel Reagent Reservoirs |
| | A1 | 25 ml Multichannel Reagent Reservoir (Insert) | INTEGRA | 4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382 | Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only |
| | A2 | 25 ml Multichannel Reagent Reservoir (Insert) | INTEGRA | 4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382 | Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only |
| B | Tube Rack | Rack for 1.5 ml microcentrifuge tubes - 1500 µl | INTEGRA | 4540 | 6x8 1.5 ml microcentrifuge tubes |
| C | Plate | 384 Well F-Bottom Low Vol. Plate - 50 µl | CORNING | 3820, 3821, 3822, 3824, 3825, 3826, 3540, 3542, 4518, 4681, 4581, 4583, 4585, 4587 | |
| D | Waste | | | | |

Method


1 Initial Volumes

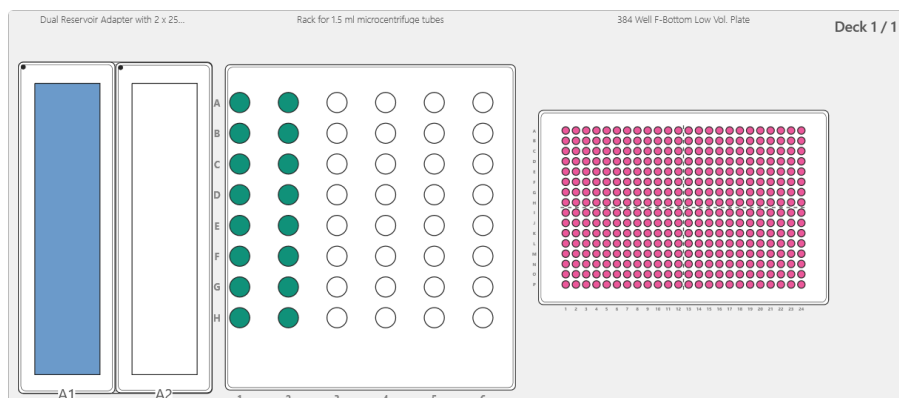

2 Repeat Dispense



Time:
2 min 1 sec
Used Tips:
8

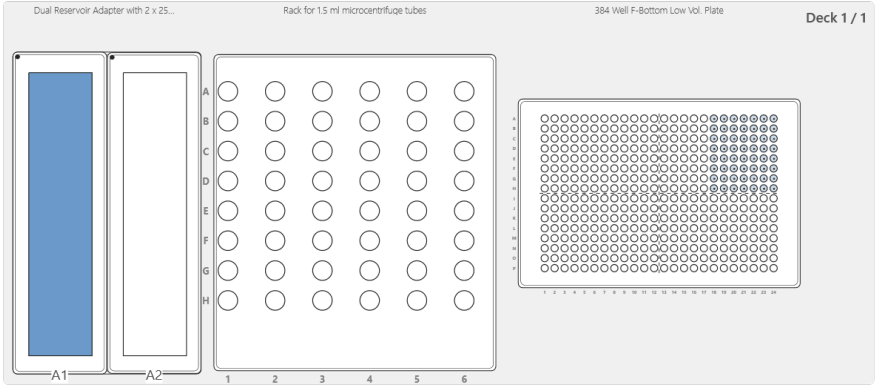

Summary individual transfers

| Step | Source | | | Target | | | Volume [µl] |
|------|---------------|----------------|-------------------|---------------|----------------|-------------------|-------------|
| | Deck Position | Well Positions | Start Height [mm] | Deck Position | Well Positions | Start Height [mm] | |
| 1 | A | 1 | 19.3 mm | C | A18-H18 | 4.3 mm | 20 |
| 2 | A | 1 | 19.3 mm | C | A19-H19 | 4.3 mm | 20 |
| 3 | A | 1 | 19.3 mm | C | A20-H20 | 4.3 mm | 20 |
| 4 | A | 1 | 19.3 mm | C | A21-H21 | 4.3 mm | 20 |
| 5 | A | 1 | 19.3 mm | C | A22-H22 | 4.3 mm | 20 |
| 6 | A | 1 | 19.3 mm | C | A23-H23 | 4.3 mm | 20 |
| 7 | A | 1 | 19.3 mm | C | A24-H24 | 4.3 mm | 20 |


Pipetting settings

| Tab | Parameter | Set value |
|------------------|--|--|
| Volumes | Volume Pre-Dispense Post-Dispense Post-Dispense Location Reuse Post-Dispense Dispense Type | Fix 5 µl 5 µl Source No Multi |
| Pipetting Speeds | Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense | 5 0 5 0 No No |
| Pipetting Height | Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset | Source: A1: Fix No A1: 2 mm Target: C: Fix No C: 1 mm |
| Tip Change | Tip Change | After step complete |
| Tip Touch | Tip Touch Type of Tip Touch Tip Touch Distance Tip Touch Height | Yes C: Side C: 1.2 mm C: 10.2 mm |

3 Repeat Dispense



Time:
1 min 2 sec
Used Tips:
8

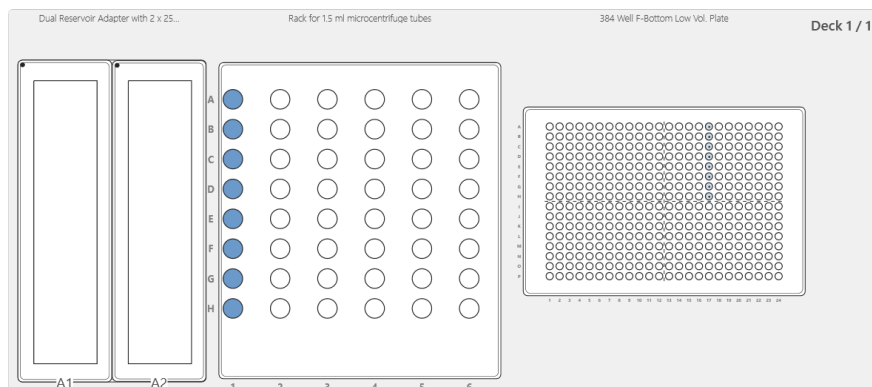

Summary individual transfers

| Step | Source | | | Target | | | Volume [µl] |
|------|---------------|----------------|-------------------|---------------|----------------|-------------------|-------------|
| | Deck Position | Well Positions | Start Height [mm] | Deck Position | Well Positions | Start Height [mm] | |
| 1 | B | A1-H1 | 14.2 mm | C | A17-H17 | 4.8 mm | 40 |

Pipetting settings

| Tab | Parameter | Set value |
|------------------|--|--|
| Volumes | Volume Pre-Dispense Post-Dispense Post-Dispense Location Reuse Post-Dispense Dispense Type | Fix 2 µl 2 µl Source No Multi |
| Pipetting Speeds | Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense | 5 0 5 0 No No |
| Pipetting Height | Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset | Source: B: Fix No B: 2 mm Target: C: Fix No C: 2 mm |
| Tip Change | Tip Change | After step complete |
| Tip Touch | Tip Touch Type of Tip Touch Tip Touch Distance Tip Touch Height | Yes C: Side C: 1.2 mm C: 10.2 mm |

Mix Summary

Mix Source

| Source | | | | |
|--------|---------------|----------------|------------------|-------------|
| Step | Deck Position | Well Positions | Pipetting Height | Volume [µl] |
| 1 | B | A1-H1 | 14.2 mm | 50 |

| Tab | Parameter | Set value |
|-----|--|--|
| Mix | Source: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel Target: Mixing | Source: Yes 6 5 1 s No Target: No |

4 Message




Pipetting settings

| Tab | Parameter | Set value |
|---------|--|------------------------------|
| Message | Message Line 1 Message Line 2 Message Line 3 | Centrifuge 1000g 5 min |

5 Serial Dilution



Time:
8 min 3 sec
Used Tips:
8

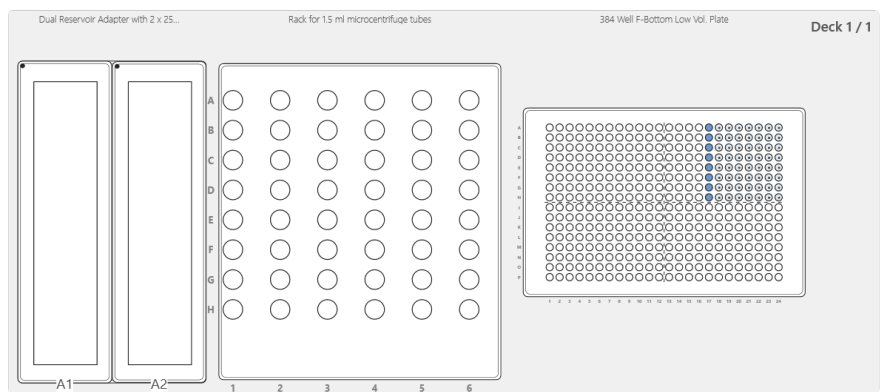

Summary individual transfers

| Step | Source | | | Target | | | Volume [µl] |
|------|---------------|----------------|-------------------|---------------|----------------|-------------------|-------------|
| | Deck Position | Well Positions | Start Height [mm] | Deck Position | Well Positions | Start Height [mm] | |
| 1 | C | A17-H17 | 4.3 mm | C | A18-H18 | 4.3 mm | 20 |
| 2 | | | | C | A19-H19 | 4.3 mm | 20 |
| 3 | | | | C | A20-H20 | 4.3 mm | 20 |
| 4 | | | | C | A21-H21 | 4.3 mm | 20 |
| 5 | | | | C | A22-H22 | 4.3 mm | 20 |
| 6 | | | | C | A23-H23 | 4.3 mm | 20 |
| 7 | | | | C | A24-H24 | 4.3 mm | 20 |

Pipetting settings

| Tab | Parameter | Set value |
|------------------|--|--|
| Volumes | Last Aspiration | Tip |
| Pipetting Speeds | Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense | 5 0 5 0 No No |
| Pipetting Height | Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset | Source: C: Fix No C: 1 mm Target: C: Fix No C: 1 mm |
| Tip Change | Tip Change | After step complete |
| Tip Touch | Tip Touch | No |

Mix Summary

Mix Target

| Target | | | | |
|--------|---------------|----------------|------------------|-------------|
| Step | Deck Position | Well Positions | Pipetting Height | Volume [µl] |
| 1 | C | - | 4.3 mm | 20 |
| 2 | C | - | 4.3 mm | 20 |
| 3 | C | - | 4.3 mm | 20 |
| 4 | C | - | 4.3 mm | 20 |
| 5 | C | - | 4.3 mm | 20 |
| 6 | C | - | 4.3 mm | 20 |
| 7 | C | - | 4.3 mm | 20 |

| Tab | Parameter | Set value |
|-----|--|---|
| Mix | Source: Mixing Target: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel | Source: No Target: Yes 12 5 4 s No |