







DLS 8 Samples Report		INTEGRAL VIALAB	
Program Name		DLS 8 Samples.iaa	
Program Name (on pipette)		DLS-8_25JAN_03	
Last Saved Date:		25. Jan 2024	
Last Save Operator:		NDziuba	
Instrument - Serial Number		0020050843	
Pipette - Serial Number		0007021615	
Tip Type (PN 6565) Lot Nr.:			
Run Operator:		NDziuba	
Run Date:		25.Jan.2024	
Run Start Time:		13:18	
Run End Time:		13:32	
Notes:		13:18:36 : Run started 13:18:39 : Repeat Dispense (Step 02) 13:19:35 : Repeat Dispense (Step 03) 13:20:32 : Wait for user input. Ready to continue protocol 13:20:32 : Message (Step 04) 13:28:12 : Run continued 13:28:13 : Serial Dilution (Step 05) 13:32:40 : Serial Dilution (Step 06) 13:32:41 : Run finished	
Signature:			
Overview Method			
<div><div><p>VOYAGER - 125µl - 8CH</p></div><div><div>1 Initial Volumes</div><div></div></div><div><div>2 Repeat Dispense</div><div></div></div><div><div>3 Repeat Dispense</div><div></div></div><div><div>4 Message</div><div></div></div><div><div>5 Serial Dilution</div><div></div></div></div>			
Total Time:		7 min 3 sec	
Total Tip Consumption:		72	
Deck Layout			

125 µl

Landscape

Dual Reservoir Adapter with 2 x 25...

Rack for 1.5 ml microcentrifuge tubes

384 Well F-Bottom Low Vol. Plate

Deck 1 / 1

A1

A2

A

25000 µl

A

B

C

D

E

F

G

H

1

2

3

4

5

6

B

1500 µl

C

50 µl

Pipette & Deck

Labware	Name	Manufacturer	Part Number
Pipette	VOYAGER 125 µl 8 channels	INTEGRA	4722
Pipette Tip	50/125 µl GripTip, Sterile, Filter, Low retention	INTEGRA	6565
Deck	3 Position Universal Deck	INTEGRA	4520

DLS 8 Samples.iaa

Page 2


25 January 2024

Deck Labware

Deck Position	Labware	Name	Manufacturer	Part Number	Description
A	COMBI System	Dual Reservoir Adapter with 2 x 25 ml Reservoirs	INTEGRA	4547	Dual Reservoir Adapter (PN 4547) with 2 x 25 ml Multichannel Reagent Reservoirs
	A1	25 ml Multichannel Reagent Reservoir (Insert)	INTEGRA	4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382	Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only
	A2	25 ml Multichannel Reagent Reservoir (Insert)	INTEGRA	4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382	Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only
B	Tube Rack	Rack for 1.5 ml microcentrifuge tubes - 1500 µl	INTEGRA	4540	6x8 1.5 ml microcentrifuge tubes
C	Plate	384 Well F-Bottom Low Vol. Plate - 50 µl	CORNING	3820, 3821, 3822, 3824, 3825, 3826, 3540, 3542, 4518, 4681, 4581, 4583, 4585, 4587	
D	Waste				

Method

1 Initial Volumes

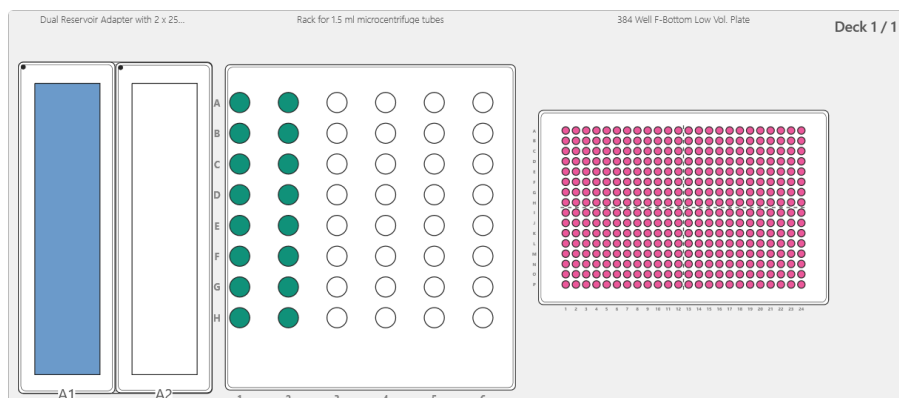


200 µl
0 µl
25000 µl

2 Repeat Dispense



Time:
1 min 12 sec
Used Tips:
8

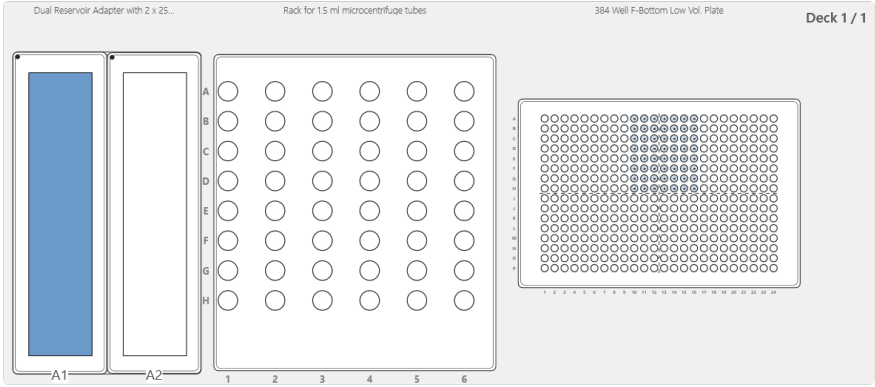

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	A	1	19.3 mm	C	A10-H10	3.8 mm	20
2	A	1	19.3 mm	C	A11-H11	3.8 mm	20
3	A	1	19.3 mm	C	A12-H12	3.8 mm	20
4	A	1	19.3 mm	C	A13-H13	3.8 mm	20
5	A	1	19.3 mm	C	A14-H14	3.8 mm	20
6	A	1	19.3 mm	C	A15-H15	3.8 mm	20
7	A	1	19.3 mm	C	A16-H16	3.8 mm	20


Pipetting settings

Tab	Parameter	Set value
Volumes	Volume Pre-Dispense Post-Dispense Post-Dispense Location Reuse Post-Dispense Dispense Type	Fix 5 µl 5 µl Source No Multi
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: A1: Fix No A1: 2 mm Target: C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch Type of Tip Touch Tip Touch Distance Tip Touch Height	Yes C: Side C: 1.2 mm C: 10.2 mm

3 Repeat Dispense



Time:
1 min 2 sec
Used Tips:
8

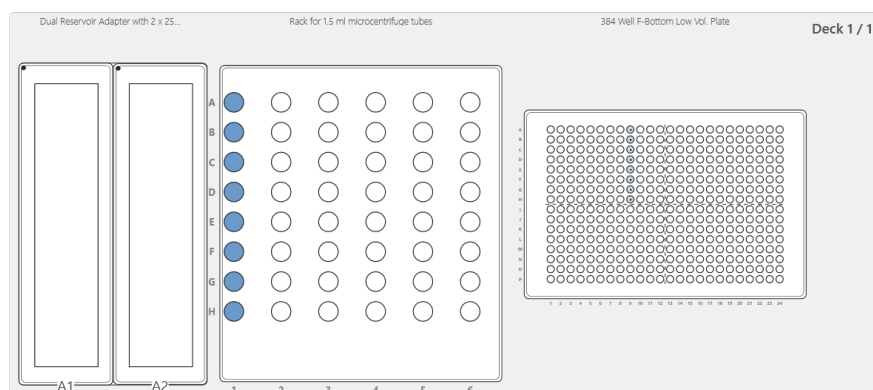

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	B	A1-H1	14.2 mm	C	A9-H9	3.8 mm	40

Pipetting settings

Tab	Parameter	Set value
Volumes	Volume Pre-Dispense Post-Dispense Post-Dispense Location Reuse Post-Dispense Dispense Type	Fix 2 µl 2 µl Source No Multi
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: B: Fix No B: 2 mm Target: C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch Type of Tip Touch Tip Touch Distance Tip Touch Height	Yes C: Side C: 1.2 mm C: 10.2 mm

Mix Summary

Mix Source

Source				
Step	Deck Position	Well Positions	Pipetting Height	Volume [µl]
1	B	A1-H1	14.2 mm	50

Tab	Parameter	Set value
Mix	Source: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel Target: Mixing	Source: Yes 6 5 1 s No Target: No

4 Message




Pipetting settings

Tab	Parameter	Set value
Message	Message Line 1 Message Line 2 Message Line 3	Centrifuge 1000g 5 min

5 Serial Dilution



Time:
4 min 43 sec
Used Tips:
56

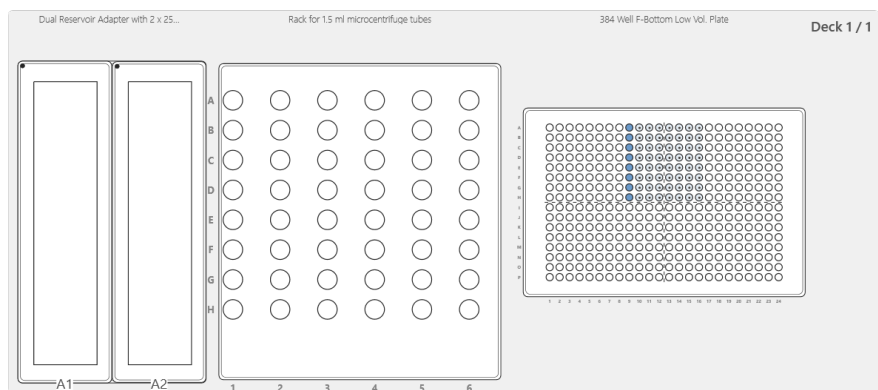

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	C	A9-H9	6.2 mm	C	A10-H10	4.3 mm	20
2				C	A11-H11	4.3 mm	20
3				C	A12-H12	4.3 mm	20
4				C	A13-H13	4.3 mm	20
5				C	A14-H14	4.3 mm	20
6				C	A15-H15	4.3 mm	20
7				C	A16-H16	4.3 mm	20

Pipetting settings

Tab	Parameter	Set value
Volumes	Last Aspiration	Tip
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: C: Fix Yes C: 1.3 mm Target: C: Fix Yes C: 1.3 mm
Tip Change	Tip Change	After each dilution
Tip Touch	Tip Touch Type of Tip Touch Tip Touch Distance Tip Touch Height	Yes C: Side C: 1.2 mm C: 10.2 mm

Mix Summary

Mix Source

Source				
Step	Deck Position	Well Positions	Start Height [mm]	Volume [μl]
1	C	-	6.4 mm	20

Mix Target

Target				
Step	Deck Position	Well Positions	Start Height [mm]	Volume [μl]
1	C	-	6.3 mm	20
2	C	-	6.3 mm	20
3	C	-	6.3 mm	20
4	C	-	6.3 mm	20
5	C	-	6.3 mm	20
6	C	-	6.3 mm	20
7	C	-	6.3 mm	20

Tab	Parameter	Set value
Mix	Source: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel Target: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel	Source: Yes 5 5 1 s Yes Target: Yes 5 5 1 s Yes