


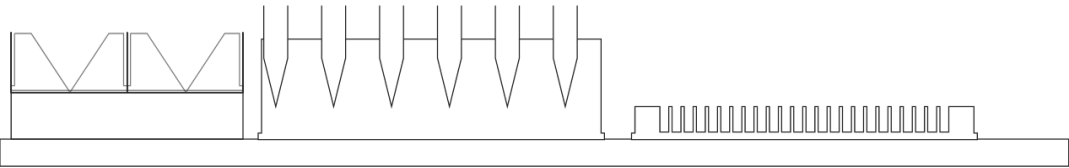
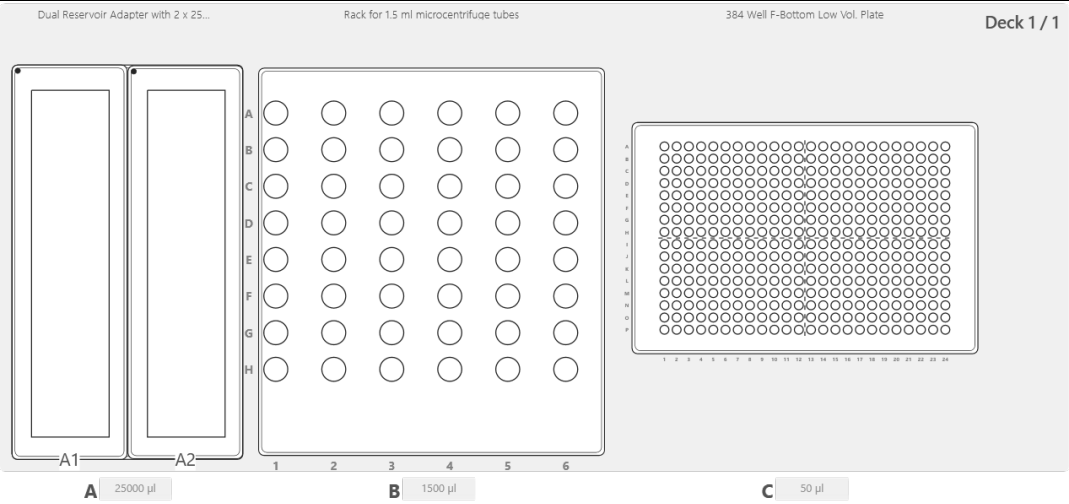
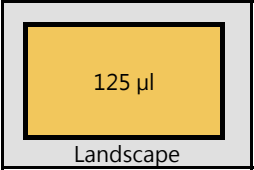


DLS 8 B 17 Dilution Series Report		INTEGRAL VIALAB	
Program Name		DLS 8 B 17 Dilution Series.iaa	
Program Name (on pipette)		DLS-8_25JAN_02	
Last Saved Date:		25. Jan 2024	
Last Save Operator:		NDziuba	
Instrument - Serial Number		0020050843	
Pipette - Serial Number		0007021615	
Tip Type (PN 6565) Lot Nr.:			
Run Operator:		NDziuba	
Run Date:		25.Jan.2024	
Run Start Time:		11:56	
Run End Time:		11:58	
Notes:		11:56:17 : Run started 11:56:21 : Serial Dilution (Step 02) 11:58:44 : Serial Dilution (Step 03) 11:58:44 : Run finished	
Signature:			
<div>Overview Method</div> <div><div><p>VOYAGER - 125µl - 8CH</p></div><div><div>1 Initial Volumes</div></div><div><div>2 Serial Dilution</div></div></div> <div><div>Total Time:</div><div>2 min 24 sec</div></div> <div><div>Total Tip Consumption:</div><div>8</div></div>			
<div>Deck Layout</div>			



Pipette & Deck


Labware	Name	Manufacturer	Part Number
Pipette	VOYAGER 125 µl 8 channels	INTEGRA	4722
Pipette Tip	50/125 µl GripTip, Sterile, Filter, Low retention	INTEGRA	6565
Deck	3 Position Universal Deck	INTEGRA	4520

Deck Labware

Deck Position	Labware	Name	Manufacturer	Part Number	Description
A	COMBI System	Dual Reservoir Adapter with 2 x 25 ml Reservoirs	INTEGRA	4547	Dual Reservoir Adapter (PN 4547) with 2 x 25 ml Multichannel Reagent Reservoirs
	A1	25 ml Multichannel Reagent Reservoir (Insert)	INTEGRA	4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382	Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only
	A2	25 ml Multichannel Reagent Reservoir (Insert)	INTEGRA	4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382	Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only
B	Tube Rack	Rack for 1.5 ml microcentrifuge tubes - 1500 µl	INTEGRA	4540	6x8 1.5 ml microcentrifuge tubes
C	Plate	384 Well F-Bottom Low Vol. Plate - 50 µl	CORNING	3820, 3821, 3822, 3824, 3825, 3826, 3540, 3542, 4518, 4681, 4581, 4583, 4585, 4587	
D	Waste				

Method


1 Initial Volumes

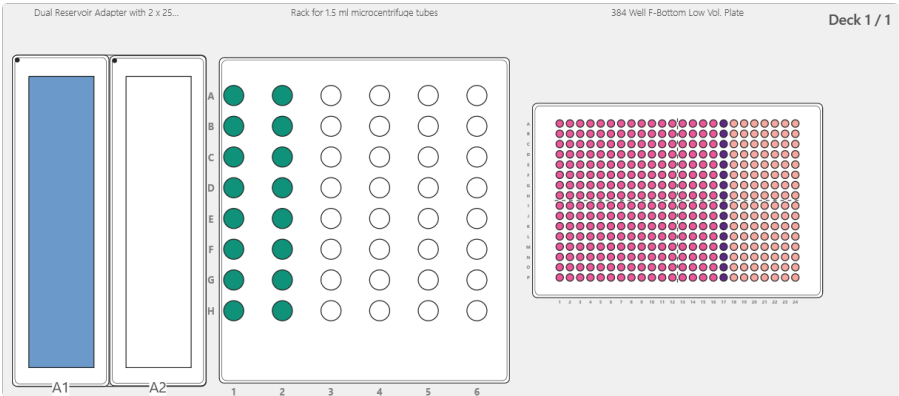

2 Serial Dilution



Time:  
2 min 19 sec  
Used Tips:  
8

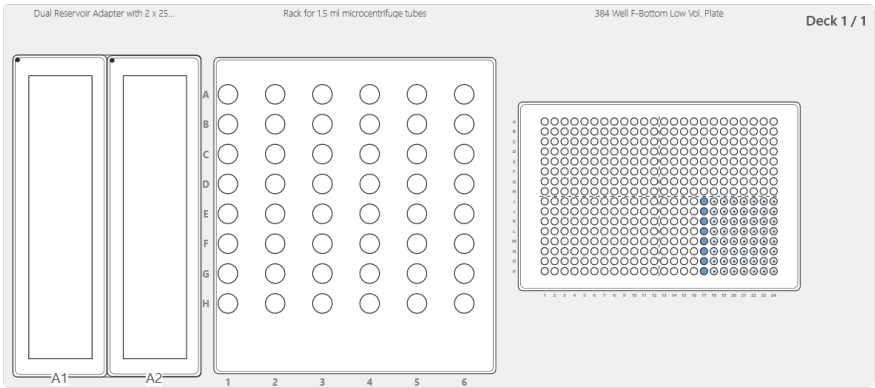

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	C	I17-P17	4.3 mm	C	I18-P18	4.3 mm	20
2				C	I19-P19	4.3 mm	20
3				C	I20-P20	4.3 mm	20
4				C	I21-P21	4.3 mm	20
5				C	I22-P22	4.3 mm	20
6				C	I23-P23	4.3 mm	20
7				C	I24-P24	4.3 mm	20

## Pipetting settings

Tab	Parameter	Set value
Volumes	Last Aspiration	Tip
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0  No No
Pipetting Height	<b>Source:</b> Heights Tip Travel Safety Bottom Offset  <b>Target:</b> Heights Tip Travel Safety Bottom Offset	<b>Source:</b> C: Fix No C: 1 mm  <b>Target:</b> C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch	No

## Mix Summary

### Mix Source

Source				
Step	Deck Position	Well Positions	Pipetting Height	Volume [µl]
1	C	-	4.3 mm	20

### Mix Target

Target				
Step	Deck Position	Well Positions	Pipetting Height	Volume [µl]
1	C	-	4.3 mm	20
2	C	-	4.3 mm	20
3	C	-	4.3 mm	20
4	C	-	4.3 mm	20
5	C	-	4.3 mm	20
6	C	-	4.3 mm	20
7	C	-	4.3 mm	20

Tab	Parameter	Set value
Mix	<b>Source:</b> Mixing Mix Cycles Mix Speed Mix Pause Tip Travel	<b>Source:</b> Yes 5 5 1 s No
	<b>Target:</b> Mixing Mix Cycles Mix Speed Mix Pause Tip Travel	<b>Target:</b> Yes 5 5 1 s No