


Program Name	DLS 8 Samples.iaa
Program Name (on pipette)	DLS-8_25JAN_01
Last Saved Date:	25. Jan 2024
Last Save Operator:	NDziuba
Instrument - Serial Number	0020050843
Pipette - Serial Number	0007021615
Tip Type (PN 6565) Lot Nr.:	
Run Operator:	NDziuba
Run Date:	25.Jan.2024
Run Start Time:	11:31
Run End Time:	11:48
Notes:	<p>11:31:59 : Run started</p> <p>11:32:02 : Repeat Dispense (Step 02)</p> <p>11:34:15 : Repeat Dispense (Step 03)</p> <p>11:35:13 : Wait for user input. Ready to continue protocol</p> <p>11:35:13 : Message (Step 04)</p> <p>11:43:45 : Run continued</p> <p>11:43:46 : Serial Dilution (Step 05)</p> <p>11:45:51 : Serial Dilution (Step 06)</p> <p>11:48:01 : Serial Dilution (Step 07)</p> <p>11:48:06 : Run pausing</p> <p>11:48:09 : Run paused</p> <p>11:48:09 : Wait for user input. Ready to continue protocol</p> <p>11:48:13 : Run paused</p> <p>11:48:20 : Run aborted (Eject Tip Go to HOME position)</p>
Signature:	

Overview Method




VOYAGER -
125µl - 8CH

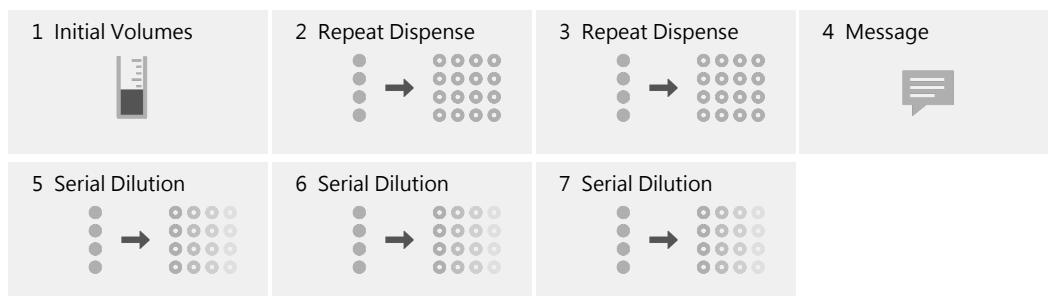

Total Time: 9 min 35 sec
Total Tip Consumption: 40

Deck Layout

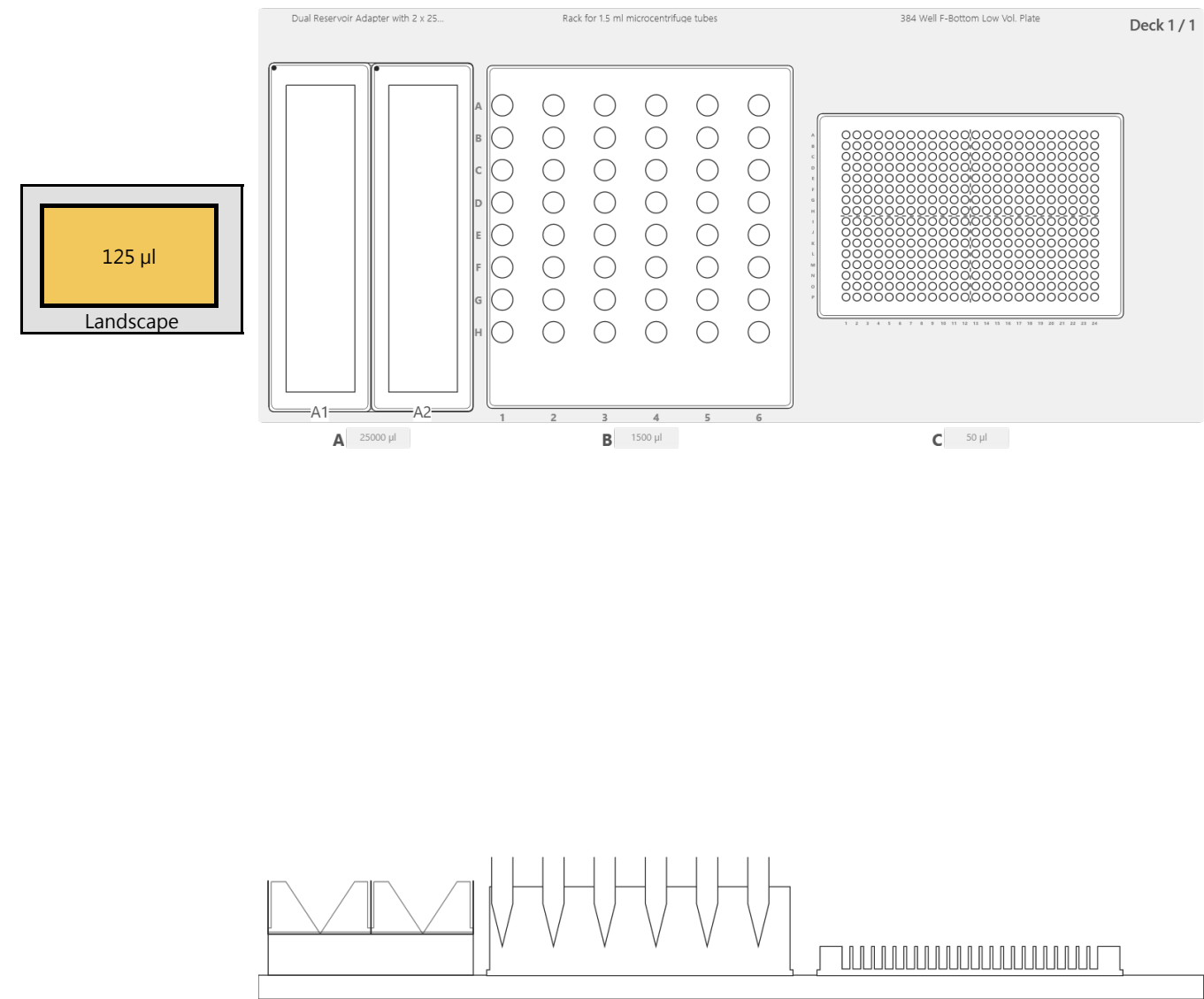

Pipette & Deck

Labware	Name	Manufacturer	Part Number
Pipette	VOYAGER 125 µl 8 channels	INTEGRA	4722
Pipette Tip	50/125 µl GripTip, Sterile, Filter, Low retention	INTEGRA	6565
Deck	3 Position Universal Deck	INTEGRA	4520

Deck Labware

Deck Position	Labware	Name	Manufacturer	Part Number	Description
A	COMBI System	Dual Reservoir Adapter with 2 x 25 ml Reservoirs	INTEGRA	4547	Dual Reservoir Adapter (PN 4547) with 2 x 25 ml Multichannel Reagent Reservoirs
	A1	25 ml Multichannel Reagent Reservoir (Insert)	INTEGRA	4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382	Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only
	A2	25 ml Multichannel Reagent Reservoir (Insert)	INTEGRA	4310, 4311, 4312, 4315, 4316, 4317, 4380, 4381, 4382	Polystyrene or Polypropylene use with Dual Reservoir Adapter (PN 4547) only
B	Tube Rack	Rack for 1.5 ml microcentrifuge tubes - 1500 µl	INTEGRA	4540	6x8 1.5 ml microcentrifuge tubes
C	Plate	384 Well F-Bottom Low Vol. Plate - 50 µl	CORNING	3820, 3821, 3822, 3824, 3825, 3826, 3540, 3542, 4518, 4681, 4581, 4583, 4585, 4587	
D	Waste				

Method

1 Initial Volumes




Deck 1 / 1

200 µl
0 µl
25000 µl

2 Repeat Dispense



Time:
2 min 45 sec
Used Tips:
8

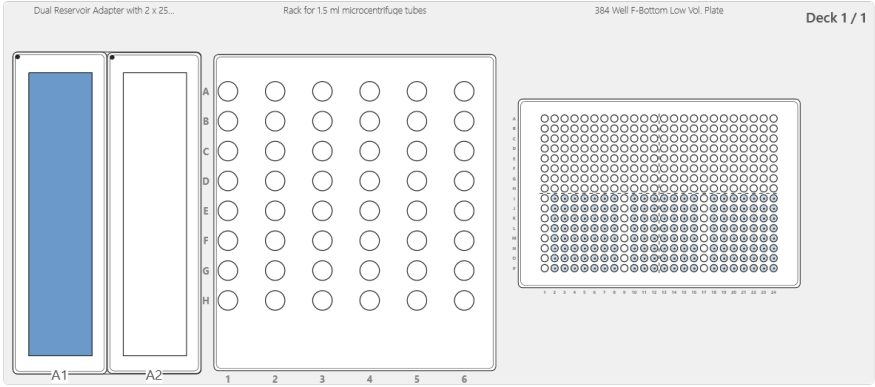

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	A	1	19.3 mm	C	I2-P2	3.8 mm	20
2	A	1	19.3 mm	C	I3-P3	3.8 mm	20
3	A	1	19.3 mm	C	I4-P4	3.8 mm	20
4	A	1	19.3 mm	C	I5-P5	3.8 mm	20
5	A	1	19.3 mm	C	I6-P6	3.8 mm	20
6	A	1	19.3 mm	C	I7-P7	3.8 mm	20
7	A	1	19.3 mm	C	I8-P8	3.8 mm	20
8	A	1	19.3 mm	C	I10-P10	3.8 mm	20
9	A	1	19.3 mm	C	I11-P11	3.8 mm	20
10	A	1	19.3 mm	C	I12-P12	3.8 mm	20
11	A	1	19.3 mm	C	I13-P13	3.8 mm	20
12	A	1	19.3 mm	C	I14-P14	3.8 mm	20
13	A	1	19.3 mm	C	I15-P15	3.8 mm	20
14	A	1	19.3 mm	C	I16-P16	3.8 mm	20
15	A	1	19.3 mm	C	I18-P18	3.8 mm	20
16	A	1	19.3 mm	C	I19-P19	3.8 mm	20
17	A	1	19.3 mm	C	I20-P20	3.8 mm	20
18	A	1	19.3 mm	C	I21-P21	3.8 mm	20
19	A	1	19.3 mm	C	I22-P22	3.8 mm	20
20	A	1	19.3 mm	C	I23-P23	3.8 mm	20
21	A	1	19.3 mm	C	I24-P24	3.8 mm	20

Pipetting settings

Tab	Parameter	Set value
Volumes	Volume Pre-Dispense Post-Dispense Post-Dispense Location Reuse Post-Dispense Dispense Type	Fix 5 µl 5 µl Source No Multi
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: A1: Fix No A1: 2 mm Target: C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch Type of Tip Touch Tip Touch Distance Tip Touch Height	Yes C: Side C: 1.2 mm C: 10.2 mm

3 Repeat Dispense




Time:

1 min 4 sec

Used Tips:

8

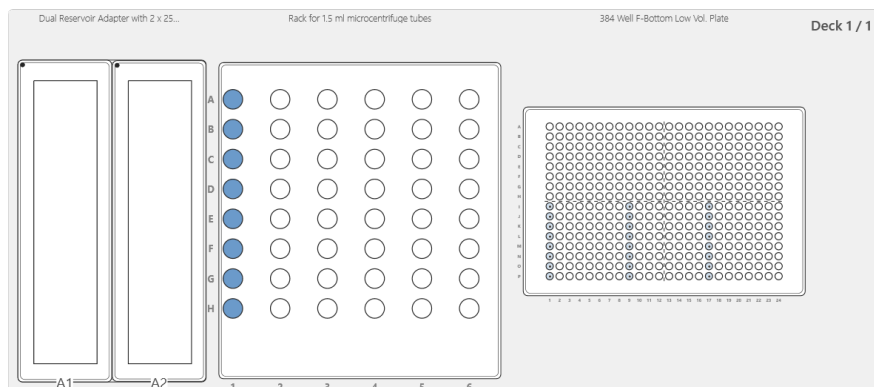

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	B	A1-H1	14.2 mm	C	I1-P1	3.8 mm	40
2	B	A1-H1	14.2 mm	C	I9-P9	3.8 mm	40
3	B	A1-H1	14.2 mm	C	I17-P17	3.8 mm	40

Pipetting settings

Tab	Parameter	Set value
Volumes	Volume Pre-Dispense Post-Dispense Post-Dispense Location Reuse Post-Dispense Dispense Type	Fix 2 µl 2 µl Source No Multi
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: B: Fix No B: 2 mm Target: C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch Type of Tip Touch Tip Touch Distance Tip Touch Height	Yes C: Side C: 1.2 mm C: 10.2 mm

Mix Summary

Mix Source

Step	Source			Volume [µl]
	Deck Position	Well Positions	Pipetting Height	
1	B	A1-H1	14.2 mm	50

Tab	Parameter	Set value
Mix	Source: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel Target: Mixing	Source: Yes 6 8 1 s No Target: No

4 Message




Pipetting settings

Tab	Parameter	Set value
Message	Message Line 1 Message Line 2 Message Line 3	Centrifuge 1000g 5 min

5 Serial Dilution



Time:
1 min 53 sec
Used Tips:
8

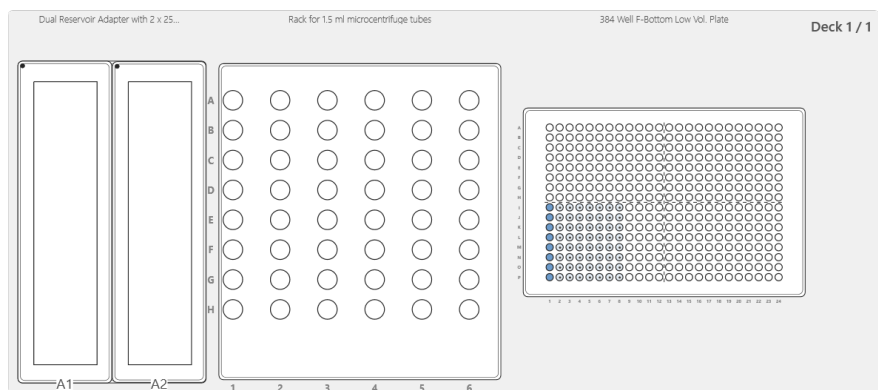

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	C	I1-P1	3.8 mm	C	I2-P2	3.8 mm	20
2				C	I3-P3	3.8 mm	20
3				C	I4-P4	3.8 mm	20
4				C	I5-P5	3.8 mm	20
5				C	I6-P6	3.8 mm	20
6				C	I7-P7	3.8 mm	20
7				C	I8-P8	3.8 mm	20

Pipetting settings

Tab	Parameter	Set value
Volumes	Last Aspiration	Tip
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: C: Fix No C: 1 mm Target: C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch	No

Mix Summary

Mix Source

Step	Source			Volume [μ l]
	Deck Position	Well Positions	Pipetting Height	
1	C	-	3.8 mm	20

Mix Target


	Target			
Step	Deck Position	Well Positions	Pipetting Height	Volume [µl]
1	C	-	3.8 mm	20
2	C	-	3.8 mm	20
3	C	-	3.8 mm	20
4	C	-	3.8 mm	20
5	C	-	3.8 mm	20
6	C	-	3.8 mm	20
7	C	-	3.8 mm	20

Tab	Parameter	Set value
Mix	Source:	Source:
	Mixing	Yes
	Mix Cycles	5
	Mix Speed	8
	Mix Pause	1 s
	Tip Travel	No
	Target:	Target:
	Mixing	Yes
	Mix Cycles	5
	Mix Speed	8
	Mix Pause	1 s
	Tip Travel	No

6 Serial Dilution



Time:
1 min 53 sec
Used Tips:
8

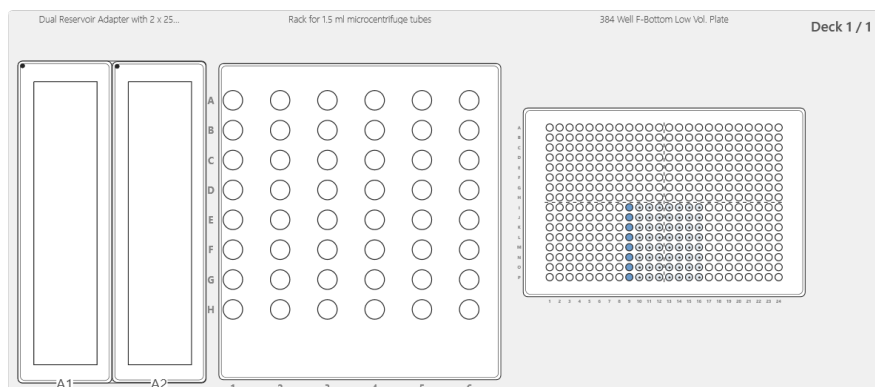

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	C	I9-P9	3.8 mm	C	I10-P10	3.8 mm	20
2				C	I11-P11	3.8 mm	20
3				C	I12-P12	3.8 mm	20
4				C	I13-P13	3.8 mm	20
5				C	I14-P14	3.8 mm	20
6				C	I15-P15	3.8 mm	20
7				C	I16-P16	3.8 mm	20

Pipetting settings

Tab	Parameter	Set value
Volumes	Last Aspiration	Tip
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: C: Fix No C: 1 mm Target: C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch	No

Mix Summary

Mix Source


Step	Source			Volume [µl]
	Deck Position	Well Positions	Pipetting Height	
1	C	-	3.8 mm	20

Mix Target

Target				
Step	Deck Position	Well Positions	Pipetting Height	Volume [µl]
1	C	-	3.8 mm	20
2	C	-	3.8 mm	20
3	C	-	3.8 mm	20
4	C	-	3.8 mm	20
5	C	-	3.8 mm	20
6	C	-	3.8 mm	20
7	C	-	3.8 mm	20

Tab	Parameter	Set value
Mix	Source: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel Target: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel	Source: Yes 5 8 1 s No Target: Yes 5 8 1 s No

7 Serial Dilution




Time:

1 min 53 sec

Used Tips:

8

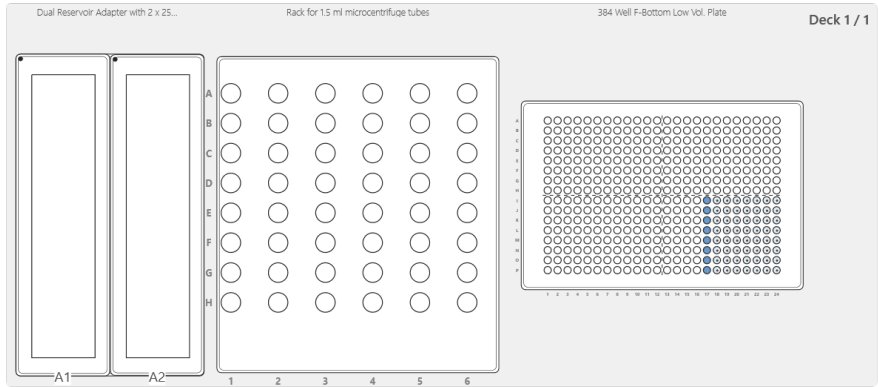

Summary individual transfers

Step	Source			Target			Volume [µl]
	Deck Position	Well Positions	Start Height [mm]	Deck Position	Well Positions	Start Height [mm]	
1	C	I17-P17	3.8 mm	C	I18-P18	3.8 mm	20
2				C	I19-P19	3.8 mm	20
3				C	I20-P20	3.8 mm	20
4				C	I21-P21	3.8 mm	20
5				C	I22-P22	3.8 mm	20
6				C	I23-P23	3.8 mm	20
7				C	I24-P24	3.8 mm	20

Pipetting settings

Tab	Parameter	Set value
Volumes	Last Aspiration	Tip
Pipetting Speeds	Aspiration Speed Aspiration Delay Dispense Speed Dispense Delay Exit Liquid Slowly Aspirate Dispense	5 0 5 0 No No
Pipetting Height	Source: Heights Tip Travel Safety Bottom Offset Target: Heights Tip Travel Safety Bottom Offset	Source: C: Fix No C: 1 mm Target: C: Fix No C: 1 mm
Tip Change	Tip Change	After step complete
Tip Touch	Tip Touch	No

Mix Summary

Mix Source

Step	Source			Volume [µl]
	Deck Position	Well Positions	Pipetting Height	
1	C	-	3.8 mm	20

Mix Target

Target				
Step	Deck Position	Well Positions	Pipetting Height	Volume [µl]
1	C	-	3.8 mm	20
2	C	-	3.8 mm	20
3	C	-	3.8 mm	20
4	C	-	3.8 mm	20
5	C	-	3.8 mm	20
6	C	-	3.8 mm	20
7	C	-	3.8 mm	20

Tab	Parameter	Set value
Mix	Source: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel Target: Mixing Mix Cycles Mix Speed Mix Pause Tip Travel	Source: Yes 5 8 1 s No Target: Yes 5 8 1 s No